


Stainless Steel Inserts



Trays, Baskets and Case Trays

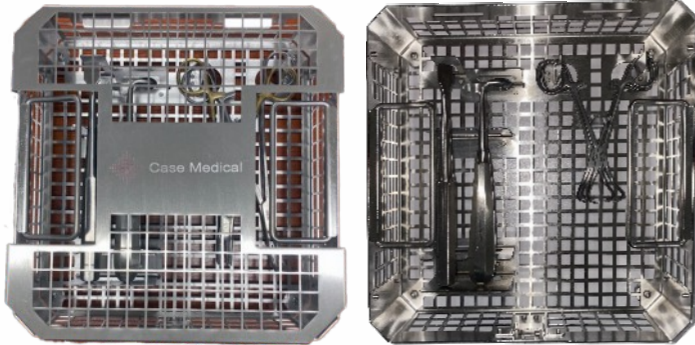

Stainless Steel Dual Use Case Tray



Stainless Steel Custom Insert with Graphics

Meditray Stainless Steel Baskets and Case Trays correspond to Case Medical's full line of SteriTite containers



Custom and Modular Inserts:

Case Medical provides inserts including trays and baskets that can be containerized or wrapped in medical grade wrappers. Aluminum inserts are commonly used because they are light in weight and highly thermoconductive especially for steam sterilization and contribute to drier loads. Stainless steel inserts are also available from Case Medical for similar use.

Dual Purpose:

Case Medical Dual Use Case Trays feature a lid, base, and locking mechanism. Use the provided screw-on feet when processed in SteriTite® and other DIN-sized sealed containers. Remove the feet to avoid tears when wrapped. Dual use case trays can be wrapped or containerized.

Custom Stainless Steel Trays and Case Trays:

Case Medical provides custom tray design utilizing advanced CAD/CAM technology, laser cutting, and laser marking for graphic layouts and labeling.

Cleaning Instructions

A pH neutral detergent or enzymatic cleaner is required to maintain the anodized seal of aluminum inserts, trays, baskets, and case trays. A mild alkaline detergent maybe used for the stainless steel inserts and case trays

